

# Wireless. Digital. Affordable.

SMART IMAGING SOLUTIONS


### The Budget-Friendly Way to Upgrade to DR Imaging.

Stepping up to full digital X-ray is a tremendous way to speed your workflow, increase productivity and support a higher standard of care. But the cost of upgrading can be a real barrier.

That's why Carestream has introduced our highly economical CARESTREAM Focus 35C Detector, powered by CARESTREAM Image Suite Software. It delivers the productivity power of wireless operation and the image quality only DR technology can provide. Best of all, the Focus Detector slides right into your existing analog equipment – protecting your current investment.

### Advantages at a Glance.

- An affordable way to upgrade to full digital imaging.
- The superb image quality of DR supports greater diagnostic confidence.
- Wireless design boosts workflow.
- Image Suite V4 Software includes a full suite of veterinary measurement tools.

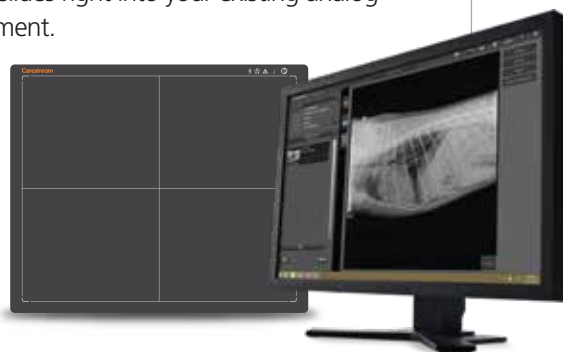



# Focus 35C Detector with Image Suite V4 Software

## Say Goodbye to the Drawbacks of Analog Imaging.

Film and Computed Radiography are both dependable imaging modalities. But the many steps of film processing can really slow you down. Plus, film means working with noxious-smelling chemicals, maintaining a darkroom and the ongoing purchase of consumables. CR imaging is consumable-free, but still requires the extra steps of scanning and digitizing each captured image, which slows down your workflow.

Our Focus Detector with Image Suite V4 Software can make all of these problems of the past. It's the ideal way to tap the power of digital – and it's budget-friendly. Additional Focus Detector advantages include:

- 35x43 cm size with Cesium Iodide scintillator for low-dose applications.
- Wireless design also helps minimize the infection risks and trip hazards of cables.
- Beam Detect Mode eliminates cable connection to the generator, alleviating OEM service concerns; image capture begins automatically when the detector senses the X-ray exposure.

## The Power Behind the Detector: Image Suite Software.

The Image Suite software platform combines advanced image processing and robust functionality for maximum efficiency. Image Suite Software:

- Provides a simple and intuitive user interface custom-designed for veterinary practices.
- Allows automated workflow with uniform protocols.
- Adds extra Cybersecurity with the Windows 10 Platform.
- Offers an optional Mini-PACS module with a broad tool set, allowing superb image management and analysis.
- Performs backups to CDs and USB drives for reliable off-line storage.
- Offers a tablet-viewing option\* for image review with patients to help you better explain your diagnosis and treatment planning.

## Image Suite Software's Advanced Reading and Measurement Tools.

Hip Dysplasia



Tibial Tuberosity Advancement



Vertebral Heart Score



\* Tablet-Viewer Software is a purchasable option and requires a separate license.

"Rx only"

© Carestream Health, Inc., 2020. CARESTREAM is a trademark of Carestream Health. CAT 2000 276 02/20

## Focus 35C Detector Specifications:

Size	35x43 cm
Weight	3.43 kg/7.56 lbs.
Battery Life	6.5 hrs
Spatial Resolution	3.3 lp/mm
Scintillator	CsI



[carestream.com](http://carestream.com)



# Carestream